

Value when  
quality counts




TriCare™

# How does TriCare™ deliver both value and comfort?

Multiple daily insulin injections can take their toll on your patients.

The TriCare™ pen needle range is manufactured with internationally approved quality standards to provide patients with a comfortable injection experience.

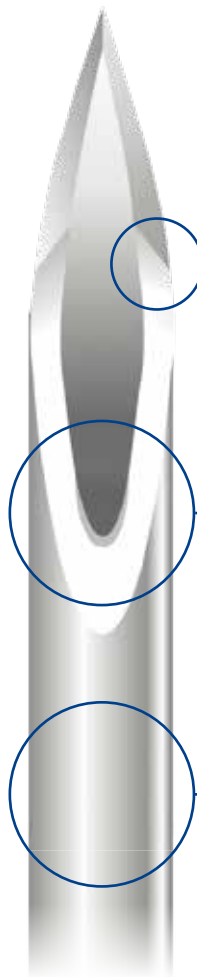


Triple sharpened bevelled edged needle tip for **clean and precise skin penetration.**


Needle body and tip are electro-polished to remove small blemishes and rough spots for a **smooth penetration.**

Silicone lubricant for **smooth gliding and less friction.**

Thin wall technology for **improved drug flow** and **reduced injection force.**

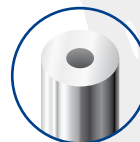

Normal-walled

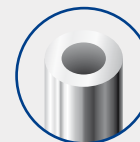

Thin-walled

- Compatible with all major brands of injection pens<sup>1</sup>.
- Non-toxic, latex free.
- Complies with ISO 11608-2.
- CE Marked.



1. Compatibility data correct as of June 2016.

## TriCare™

For more information please visit us at [owenmumford.com](http://owenmumford.com)